



OrangeDot Face Shield Version 3B Quick Start Guide 4.10.2020

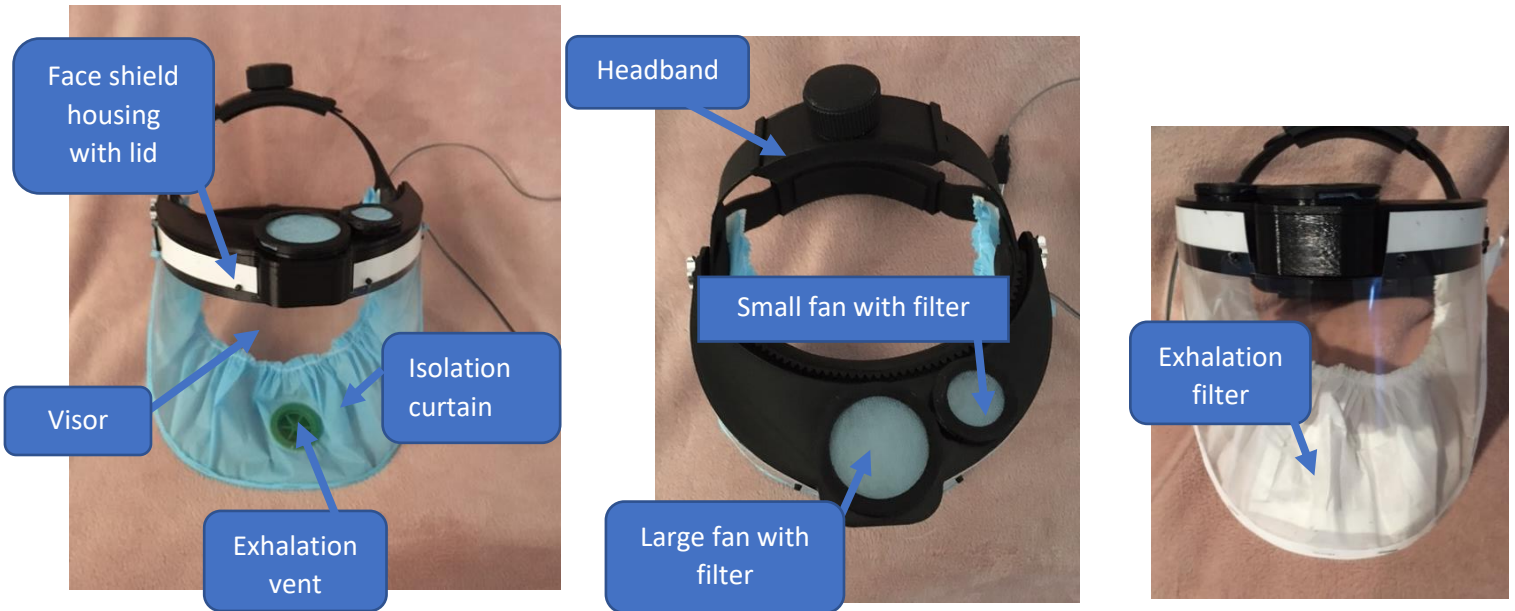

- Please undertake the Safety Test before using this device; see instructions that came with your device or visit www.orangedot.com.au.

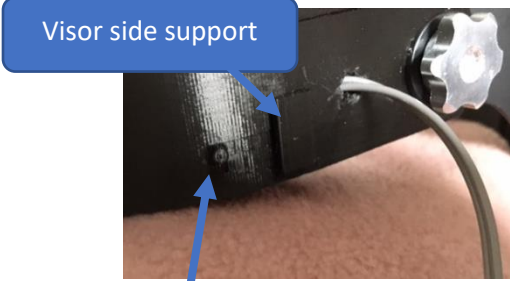

- Carefully attach the visor with isolation curtain to the housing. The edge of the visor should be tucked into the visor side support and you will need to use some force to push the visor onto the four retaining posts.
- Gently pull the free ends of the curtain and attach them to the Velcro type retaining strips on the inside of the headband. This will help to maintain a seal around the face.

Visor retaining post



Retaining strips

Headband



Curtain attached to the inside of headband

Isolation curtain



Aluminium wheel

Power cable

- Attach the USB power cord to a power bank. The two fans will start immediately. The power bank can be kept in your pocket. A 2,200 mAh power bank should last 6 hours.
- Gently insert your head into the face shield chamber.
- Adjust the vertical and horizontal headband black control wheels, so that the shield fits comfortably on your head.

elastic band backwards, so that it forms a seal around your face. Remember that this seal will not be airtight.

- The Aluminium wheel on the side of the faced shield is non-functional in this version.
- Please do not lie the face shield on the visor- this will cause creases in the visor and the shield may fall under its own weight
- Removing the power supply to the shield or the energy depletion of the battery will result in a warning flash of the LED light and an alarm noise. Please remove the device and check your connection or replace the battery.



Inferior surface of housing

LED power off warning light



Vertical control wheel

Visor

Horizontal control



- Any questions? Please contact us admin@orangedot.com.au
- www.orangedot.com.au